
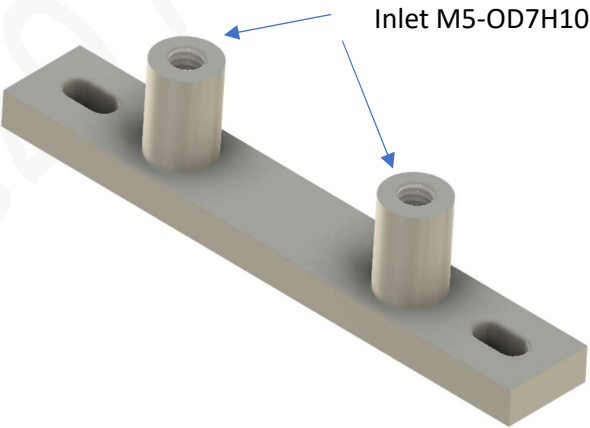



Disposable Micro Channel Air Spray Head

Drawing:

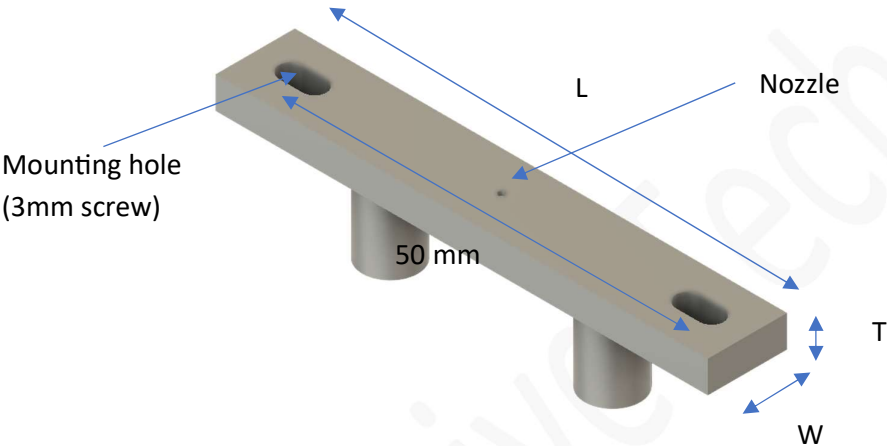

Connections:

- **Hollow screw**

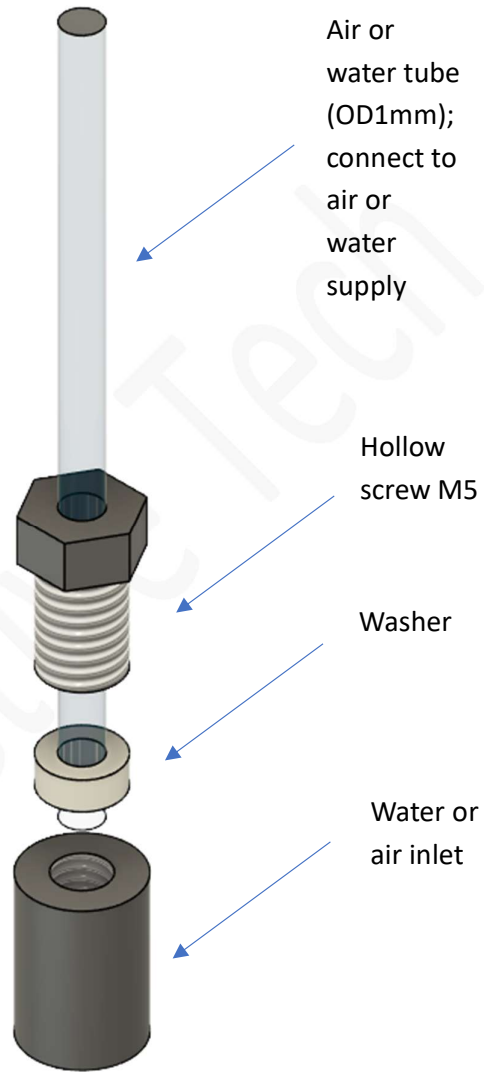

- **Luer connector / adapter / converter, etc:** One end to the M5 inlet, the other end to the tubing

Operations:

- 1) Connect the air and water tubing respectively
- 2) Turn on the water and air supply respectively
- 3) Tuning the flowrate of water and air respectively

Specifications:

Model	Materials	Nozzle	Size	Air flow rate	Water flow rate	Self-priming
MCS-621004	Acrylic	0.1mm	L62-W10-T4 In7x10M5	>3.5 L/min	Up to 5mL/min	better
MCS-621002	Acrylic	0.1mm	L62-W10-T2 In7x10M5	>2.5 L/min	Up to 5mL/min	No/weak